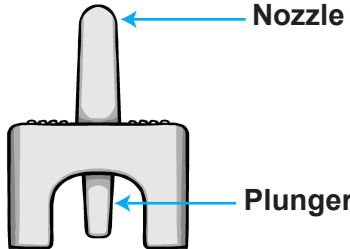
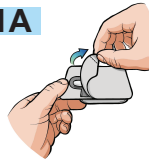

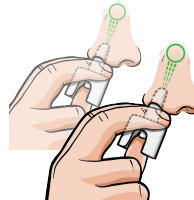
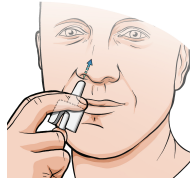
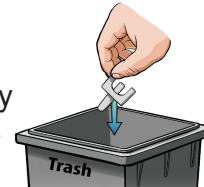
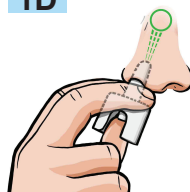




INSTRUCTIONS FOR USE NEFFY® (ne'fee) epinephrine nasal spray (Single-dose) For Use in the Nose Only	How to use the neffy nasal spray device		Get emergency medical help if needed
<p><b>Neffy nasal spray device is for allergic emergencies (anaphylaxis).</b></p> <p>Read this Instructions for Use on how to use neffy the right way. Anyone who may need to give a dose of neffy should also understand how to use the device.</p> <p>Ask a healthcare provider if you have any questions about how to use neffy.</p> <p>Always carry 2 neffy nasal spray devices with you. You may need a second dose of neffy.</p> <p><b>Your neffy nasal spray device:</b></p> <div></div> <p><b>Important Information:</b></p> <ul style="list-style-type: none"><li>•<b>Use only in the nose.</b> Avoid spraying in the eyes or mouth.</li><li>•<b>Do not</b> test or prime (pre-spray) the device. Each neffy can only be sprayed 1 time.</li><li>•Each neffy nasal spray device has only 1 dose of medicine and is for single-use only.</li><li>•Check your neffy periodically to be sure the expiration date has not passed. Replace your neffy before the expiration date has passed.</li><li>•<b>Get emergency medical help for further treatment of the allergic emergency (anaphylaxis), if needed, after using neffy.</b></li></ul>	<p><b>1A</b></p> 	<p><b>Remove neffy from packaging (see Figure 1A).</b></p> <p>Pull open the packaging to remove the neffy nasal spray device.</p>	<p><b>Get emergency medical help for further treatment of the allergic emergency (anaphylaxis), if needed, after using neffy.</b> Tell them that a neffy nasal spray has been used.</p> <p>Before you receive neffy, your healthcare provider should talk to you about when to get emergency medical help.</p>
	<p><b>1B</b></p> 	<p><b>Hold device as shown (see Figure 1B).</b></p> <p>Hold the device with your thumb on the bottom of the plunger and a finger on either side of the nozzle.</p> <ul style="list-style-type: none"><li>•<b>Do not pull or push on the plunger.</b></li><li>•<b>Do not test or prime (pre-spray).</b></li></ul> <p><b>Each device has only 1 spray.</b></p>	<p><b>After 5 minutes, repeat dose if needed</b></p> <div></div> <p>Check for symptoms after the first dose of neffy.</p> <p>If symptoms get worse or come back, give a second dose of neffy in the <b>same nostril</b> starting 5 minutes after the first dose.</p>
	<p><b>1C</b></p> 	<p><b>Insert the nozzle into a nostril until your fingers touch your nose (see Figure 1C).</b></p> <p><b>Keep the nozzle straight</b> into the nose pointed toward your forehead.</p> <p><b>Do not</b> point (angle) the nozzle to the nasal septum (wall between your 2 nostrils) or outer wall of the nose.</p>	<p><b>Throwing away the used neffy nasal spray device:</b></p> <ul style="list-style-type: none"><li>•Throw away (dispose of) the used neffy nasal spray device in household trash.</li><li>•Do not recycle.</li></ul> 
	<p><b>1D</b></p> 	<p><b>Press plunger up firmly until it snaps up and sprays liquid into the nostril (see Figure 1D).</b></p> <p><b>Do not</b> sniff during or after the dose is given. If any liquid drips out of the nose, you may need to give a second dose of neffy after checking for symptoms.</p>	<p><b>How should I store the neffy nasal spray device?</b></p> <ul style="list-style-type: none"><li>•Store neffy at room temperature between 68°F to 77°F (20°C to 25°C).</li><li>•<b>Do not freeze.</b> If neffy freezes, the device will not spray.</li><li>•Storage of neffy at high temperatures up to 122°F (50° C) is allowed for a few days.</li></ul>
	<p><b>⚠ Do not angle nozzle</b></p> <p><b>Do not</b> point (angle) the nozzle to the septum wall or outer wall of the nose.</p> <div></div>		



**neffy**<sup>®</sup>  
(epinephrine nasal spray)